



TWIST

UNI-DIRECTIONAL
DEFLECTABLE
GUIDING SHEATH


DESTINO™
TWIST

A New Twist On Access Technology



Maximizing Sheath Control

- 180° deflectable distal sheath tip section allows precise steering and navigation through the vascular system
- Clear view of the distal marker band under fluoroscopy for precise positioning

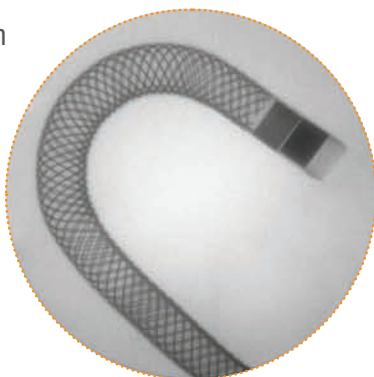


PTFE Liner

Flexsteer™ reinforced braid

Outer wall extrusion

Hydrophobic coated sheath

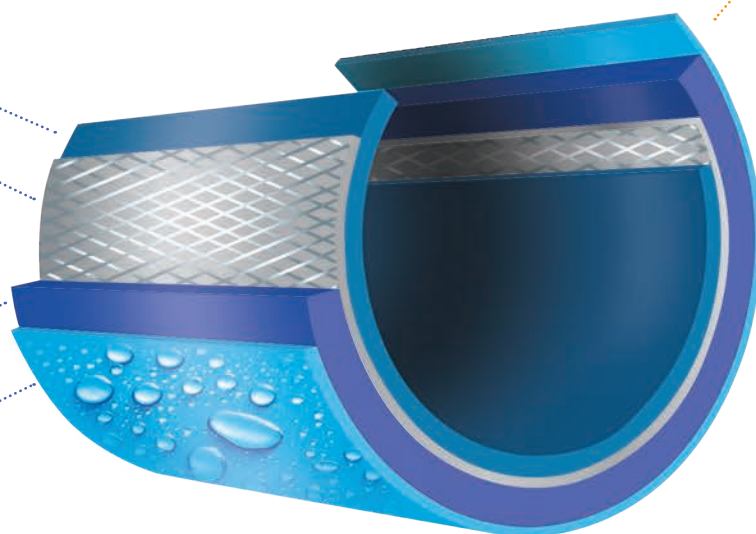


- Flexsteer™ reinforced braid technology provides exceptional torque control and kink resistance, while maintaining optimal flexibility even in the most tortuous anatomy

Minimizing Sheath Exchanges

YOUR ONE SHEATH SOLUTION

- The Destino Twist deflectable tip minimizes your need for multiple fixed curves and catheter exchanges

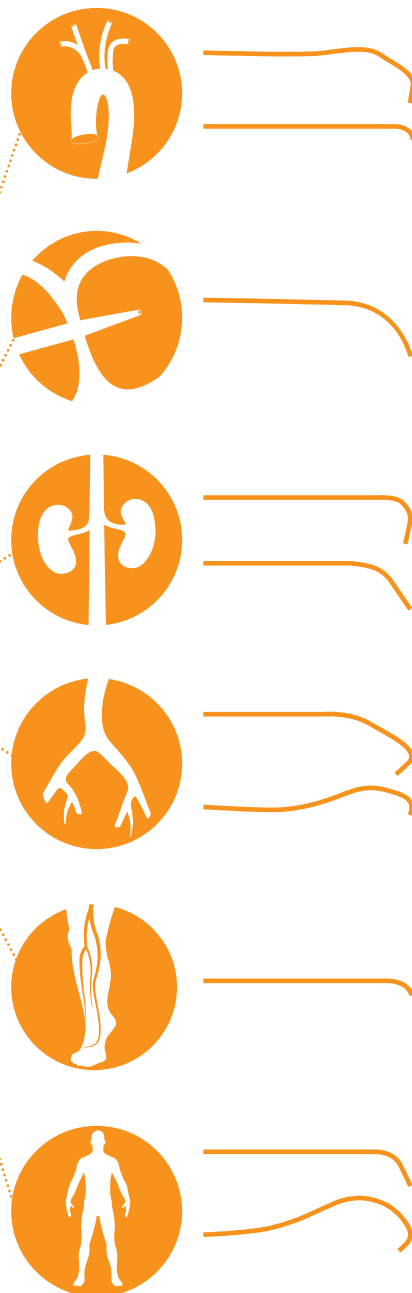


- SureSeal™ hemostatic valve technology minimizes air embolism and blood loss




Maximizing Your Options

Diameter Options:

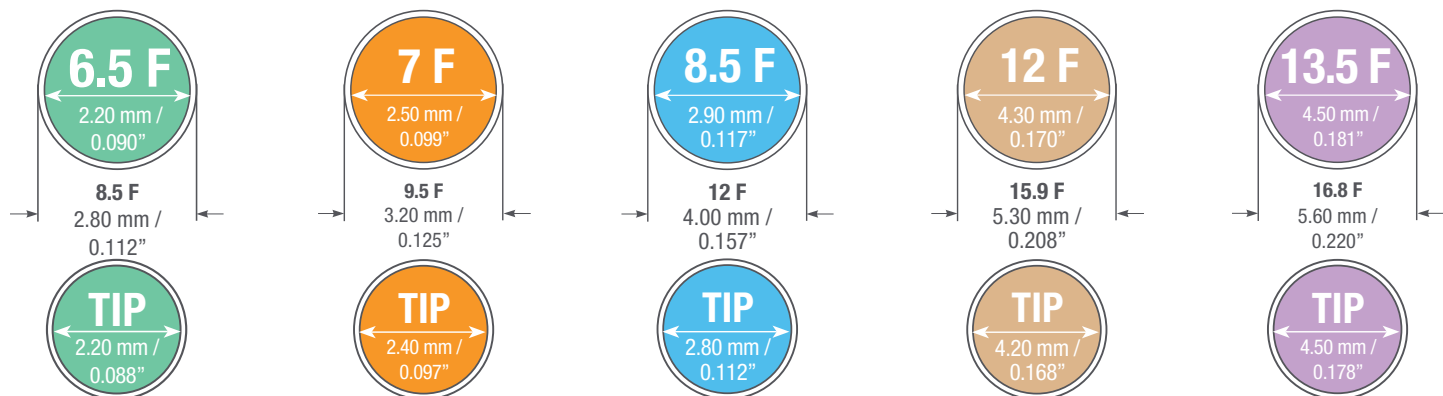

Curve Options:




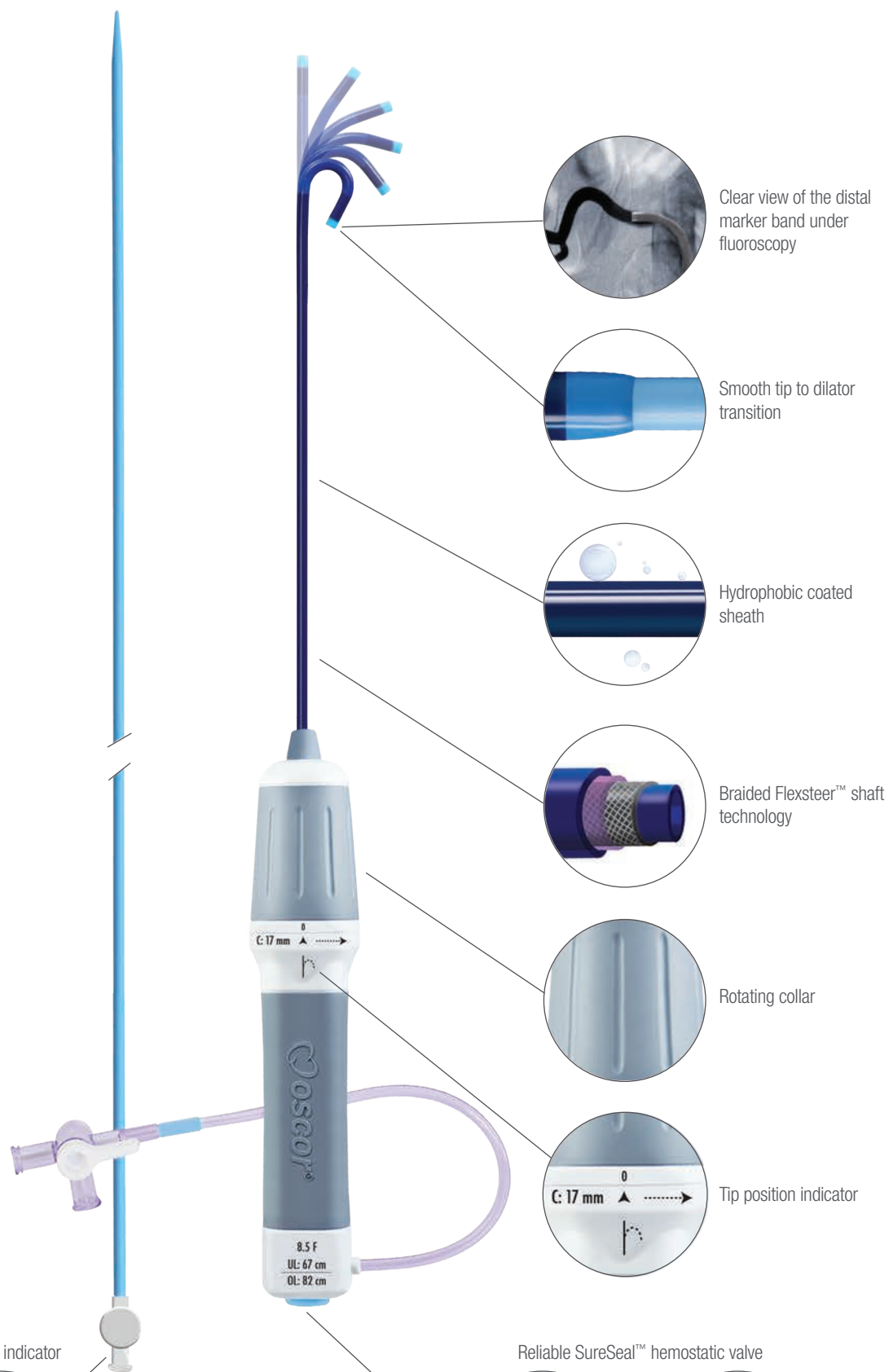

Dilator with French size and guidewire indicator



Reliable SureSeal™ hemostatic valve

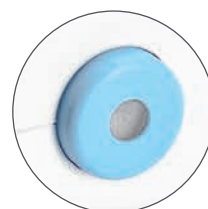





CAROTID ACCESS

A New Twist on Stability.

The Destino Twist simplifies access to challenging arch anatomy. The deflectable tip gives the user a stable platform and minimizes device exchanges.

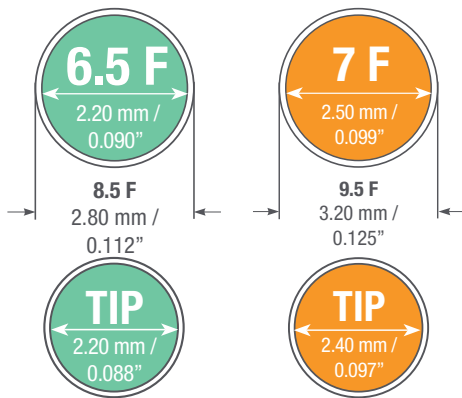


Model # Set <i>(One per box)</i>		Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend
STANDARD					
DST0659009	●	6.5 F	90 cm	108 cm	9 mm
DST0709009	●	7 F	90 cm	108 cm	9 mm

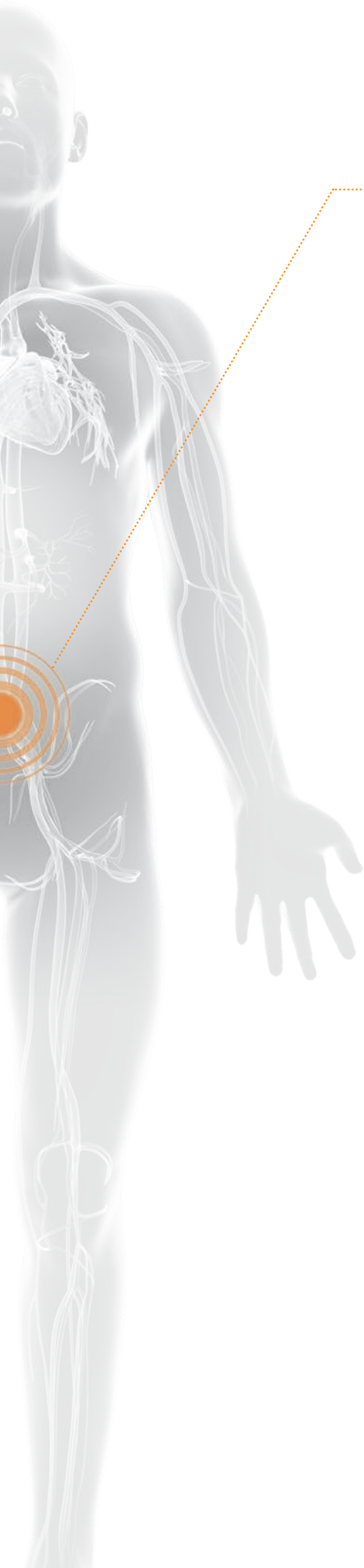

ILIAC/SFA ACCESS

A New Twist on Advanced Support.

The Destino Twist uni-directional sheath features a variety of curves, atraumatic tip, and hydrophobic coating which improves trackability. The Twist's ability to deflect, anchors and stabilizes the sheath across the bifurcation enhancing contralateral access. This also facilitates ipsilateral access of the internal iliac.

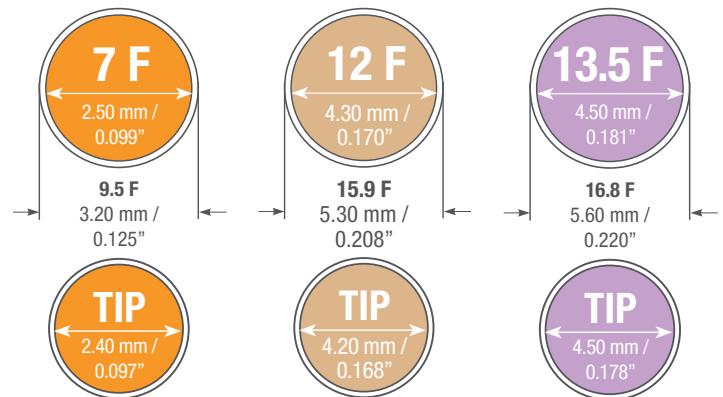

Model # Set <i>(One per box)</i>	Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend
STANDARD				
DST0704509	● 7 F	45 cm	63 cm	9 mm
DST0704517	● 7 F	45 cm	63 cm	17 mm
DST0704522	● 7 F	45 cm	63 cm	22 mm
DST0704530	● 7 F	45 cm	63 cm	30 mm
DST1204522	● 12 F	45 cm	63 cm	22 mm
DST1204530	● 12 F	45 cm	63 cm	30 mm
DST1357139	● 13.5 F	71 cm	89 cm	39 mm




INFRAINGUINAL/ INFRAPOPLITEAL ACCESS

A New Twist on Deflection Control.

The Destino Twist offers stable platform for access below the knee.

Used with a coaxial selective catheter, it simplifies access while minimizing device exchanges.

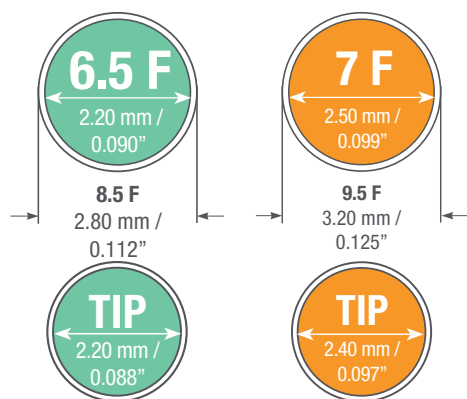


Model # Set <i>(One per box)</i>	Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend
STANDARD				
DST0654509	6.5 F	45 cm	63 cm	9 mm
DST0654517	6.5 F	45 cm	63 cm	17 mm
DST0704509	7 F	45 cm	63 cm	9 mm
DST0704517	7 F	45 cm	63 cm	17 mm



RENAL/CELIAC/SMA ACCESS

A New Twist on Accessing Challenging Anatomy.

The Destino Twist's deflectable tip simplifies access to difficult patient anatomies, which is especially useful when accessing arteries with inferior take offs. This enables the physician to access from the groin instead of using a brachial approach.

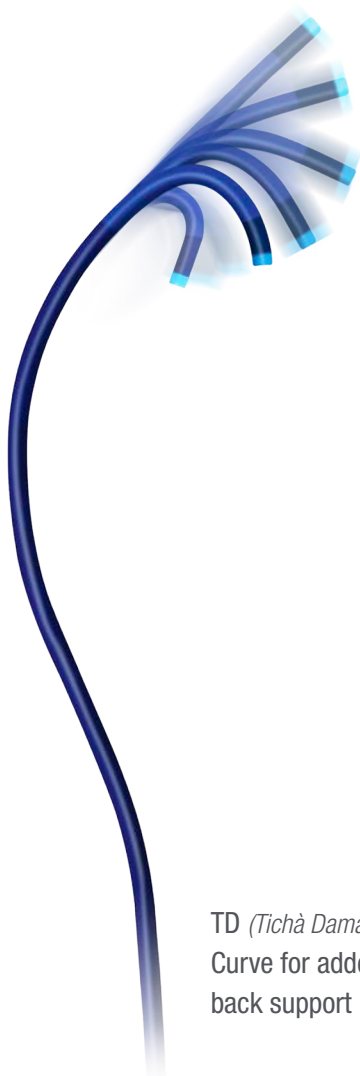

TD (*Tichà Damascelli*)
Curve for added
back support


Model # Set <i>(One per box)</i>		Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend
STANDARD					
DST0654509	●	6.5 F	45 cm	63 cm	9 mm
DST0654517	●	6.5 F	45 cm	63 cm	17 mm
DST0654522	●	6.5 F	45 cm	63 cm	22 mm
DST0654530	●	6.5 F	45 cm	63 cm	30 mm
DST0655509	●	6.5 F	55 cm	73 cm	9 mm
DST0655517	●	6.5 F	55 cm	63 cm	17 mm
DST0655522	●	6.5 F	55 cm	63 cm	22 mm
DST0655530	●	6.5 F	55 cm	63 cm	30 mm
DST0704509	●	7 F	45 cm	63 cm	9 mm
DST0704517	●	7 F	45 cm	63 cm	17 mm
DST0704522	●	7 F	45 cm	63 cm	22 mm
DST0704530	●	7 F	45 cm	63 cm	30 mm
DST0705509	●	7 F	55 cm	73 cm	9 mm
DST0705517	●	7 F	55 cm	73 cm	17 mm
DST0705522	●	7 F	55 cm	73 cm	22 mm
DST0705530	●	7 F	55 cm	73 cm	30 mm
TD					
DST0655509TD	●	6.5 F	55 cm	73 cm	9 mm
DST0705509TD	●	7 F	55 cm	73 cm	9 mm




TRANSSEPTAL ACCESS

A New Twist on Precise Advancement.

The Destino Twist features a proprietary braided Flexsteer™ hydrophobic coated sheath and the latest SureSeal™ technology, allowing for the percutaneous introduction of various types of cardiovascular catheters to all chambers of the heart, as well as transseptal needles.

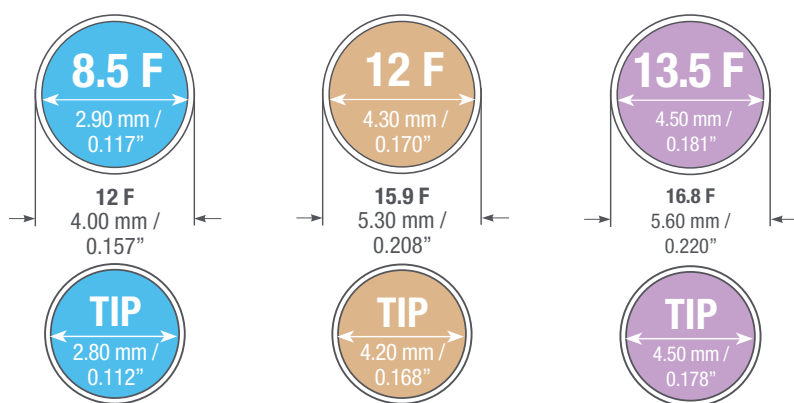


Model # Set <i>(One per box)</i>	Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend	Compatible Transseptal Needle Length
STANDARD					
DST0857317	8.5 F	73 cm	90 cm	17 mm	98 cm
DST0857322	8.5 F	73 cm	90 cm	22 mm	98 cm
DST1207122	12 F	71 cm	89 cm	22 mm	98 cm
DST1357139	13.5 F	71 cm	89 cm	39 mm	98 cm

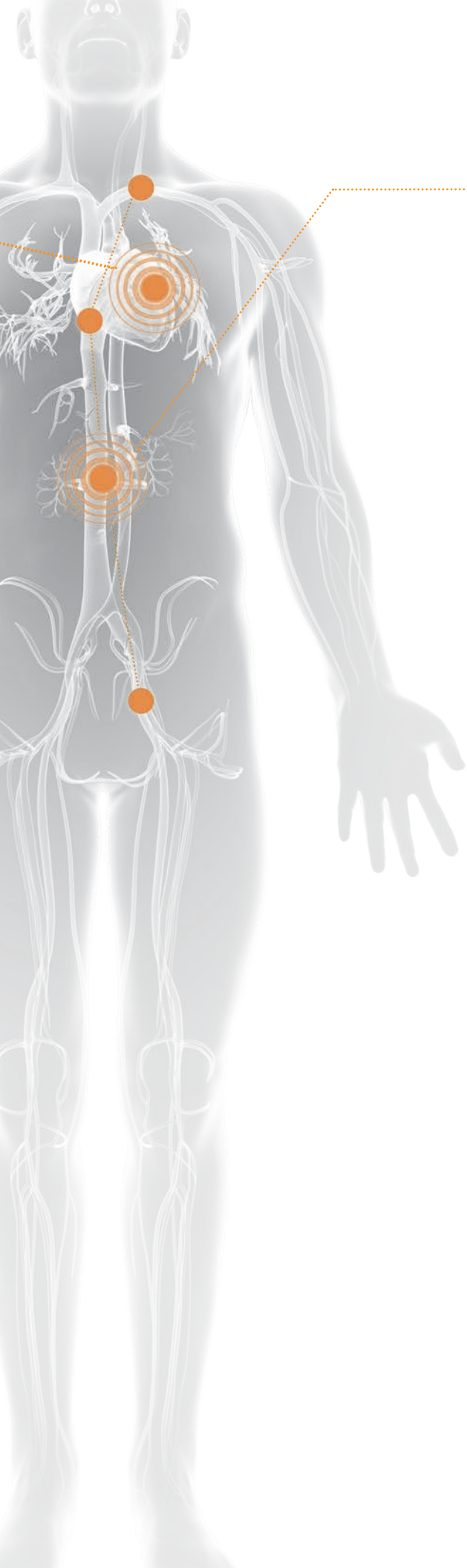

MULTIPURPOSE ACCESS

A New Twist on Adaptable Deflection.


The Destino Twist diverse product offering allows for optimal compatibility with patient anatomy adapting to all your needs.

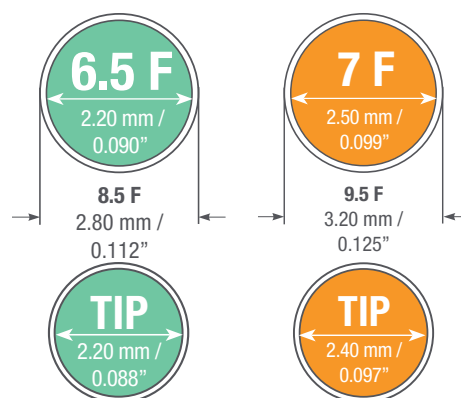

Model # Set <i>(One per box)</i>		Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend
STANDARD					
DST0655507D2*	●	6.5 F	55 cm	73 cm	7 mm
DST0655509	●	6.5 F	55 cm	73 cm	9 mm
DST0659009	●	6.5 F	90 cm	108 cm	9 mm
DST0705509	●	7 F	55 cm	73 cm	9 mm
DST0705517	●	7 F	55 cm	73 cm	17 mm
DST0709009	●	7 F	90 cm	108 cm	9 mm
DST0855517	●	8.5 F	55 cm	73 cm	17 mm
DST1207122	●	12 F	71 cm	89 cm	22 mm
DST1357139	●	13.5 F	71 cm	89 cm	39 mm

* Contains a second soft dilator with 80 cm usable length.

ORDERING INFORMATION

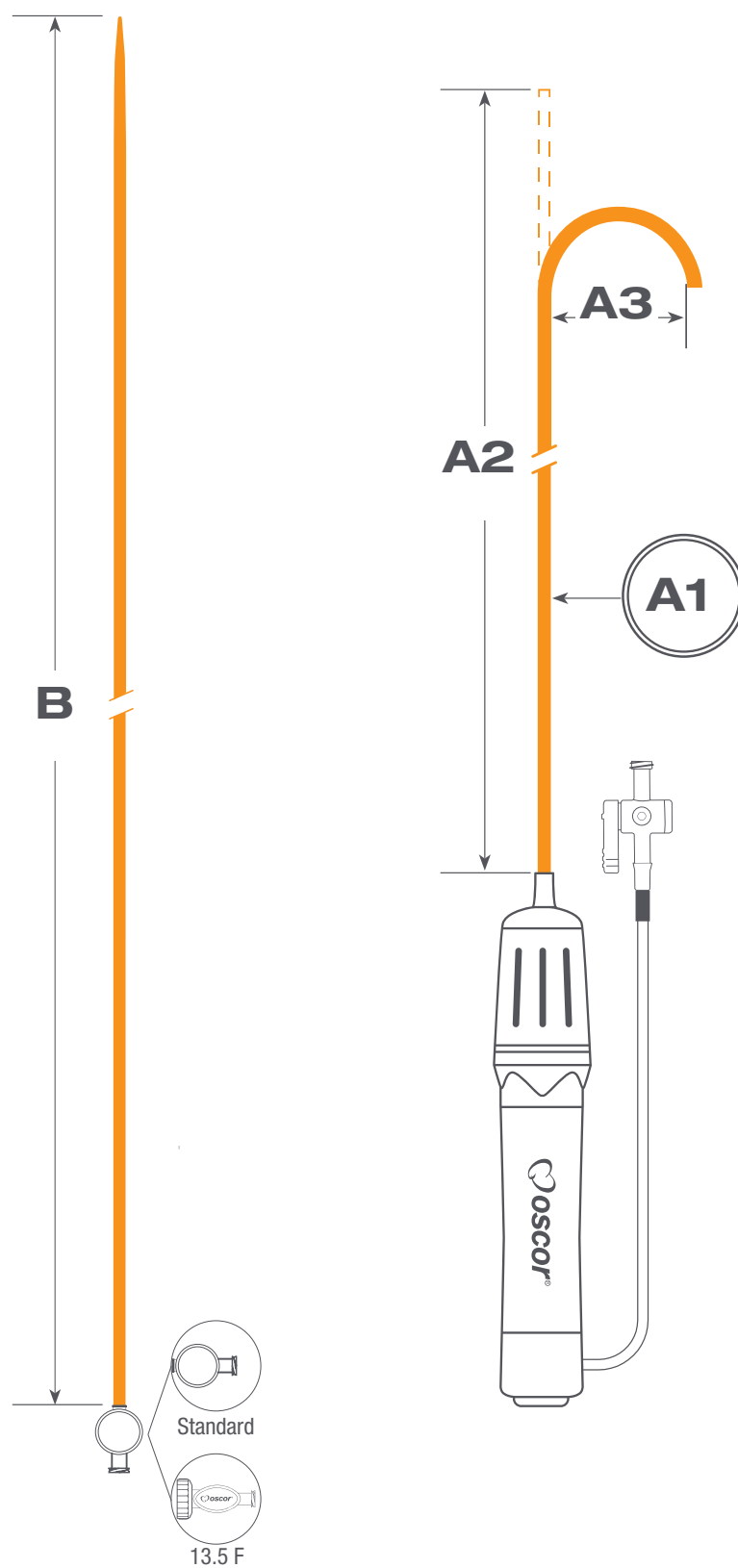



CAROTID ACCESS

Model # Set (One per box)	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD				
DST0659009	● 6.5 F	90 cm	108 cm	9 mm
DST0709009	● 7 F	90 cm	108 cm	9 mm



ILIAC / SFA ACCESS

Model # Set (One per box)	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD				
DST0704509	● 7 F	45 cm	63 cm	9 mm
DST0704517	● 7 F	45 cm	63 cm	17 mm
DST0704522	● 7 F	45 cm	63 cm	22 mm
DST0704530	● 7 F	45 cm	63 cm	30 mm
DST1204522	● 12 F	45 cm	63 cm	22 mm
DST1204530	● 12 F	45 cm	63 cm	30 mm
DST1357139	● 13.5 F	71 cm	89 cm	39 mm



INFRAINGUINAL / INFRAPOPLITEAL ACCESS

Model # Set (One per box)	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD				
DST0654509	● 6.5 F	45 cm	63 cm	9 mm
DST0654517	● 6.5 F	45 cm	63 cm	17 mm
DST0704509	● 7 F	45 cm	63 cm	9 mm
DST0704517	● 7 F	45 cm	63 cm	17 mm



RENAL / CELIAC / SMA ACCESS

Model # Set <i>(One per box)</i>	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD				
DST0654509	● 6.5 F	45 cm	63 cm	9 mm
DST0654517	● 6.5 F	45 cm	63 cm	17 mm
DST0654522	● 6.5 F	45 cm	63 cm	22 mm
DST0654530	● 6.5 F	45 cm	63 cm	30 mm
DST0655509	● 6.5 F	55 cm	73 cm	9 mm
DST0655517	● 6.5 F	55 cm	63 cm	17 mm
DST0655522	● 6.5 F	55 cm	63 cm	22 mm
DST0655530	● 6.5 F	55 cm	63 cm	30 mm
DST0704509	● 7 F	45 cm	63 cm	9 mm
DST0704517	● 7 F	45 cm	63 cm	17 mm
DST0704522	● 7 F	45 cm	63 cm	22 mm
DST0704530	● 7 F	45 cm	63 cm	30 mm
DST0705509	● 7 F	55 cm	73 cm	9 mm
DST0705517	● 7 F	55 cm	73 cm	17 mm
DST0705522	● 7 F	55 cm	73 cm	22 mm
DST0705530	● 7 F	55 cm	73 cm	30 mm
TD				
DST0655509TD	● 6.5 F	55 cm	73 cm	9 mm
DST0705509TD	● 7 F	55 cm	73 cm	9 mm



TRANSSEPTAL ACCESS

Model # Set (One per box)		A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend	Compatible Transseptal Needle Length
STANDARD						
DST0857317	●	8.5 F	73 cm	90 cm	17 mm	98 cm
DST0857322	●	8.5 F	73 cm	90 cm	22 mm	98 cm
DST1207122	●	12 F	71 cm	89 cm	22 mm	98 cm
DST1357139	●	13.5 F	71 cm	89 cm	39 mm	98 cm



MULTI PURPOSE ACCESS

Model # Set (One per box)		A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD					
DST0655507D2*	●	6.5 F	55 cm	73 cm	7 mm
DST0655509	●	6.5 F	55 cm	73 cm	9 mm
DST0659009	●	6.5 F	90 cm	108 cm	9 mm
DST0705509	●	7 F	55 cm	73 cm	9 mm
DST0705517	●	7 F	55 cm	73 cm	17 mm
DST0709009	●	7 F	90 cm	108 cm	9 mm
DST0855517	●	8.5 F	55 cm	73 cm	17 mm
DST1207122	●	12 F	71 cm	89 cm	22 mm
DST1357139	●	13.5 F	71 cm	89 cm	39 mm

* Contains a second soft dilator with 80 cm usable length.

Guiding sheaths are supplied sterile. For single use only. Available bulk upon request.
Standard lengths are noted above. Other lengths available upon request.



OSCOR INC. | WORLD HEADQUARTERS

Toll Free: 1-800-726-7267

P: +1 727-937-2511

OSCOR EUROPE GMBH | GERMANY DIRECT SALES

Telefon kostenlos: 0800 00-67267

P: +49 (0) 211 586 786-00

CONTACT US

www.oscor.com

sales@oscor.com



DST100-04062016-00